
Don't Have The Money Nor Time to Wait? Mouse Hb α , β with Human sickle Hb

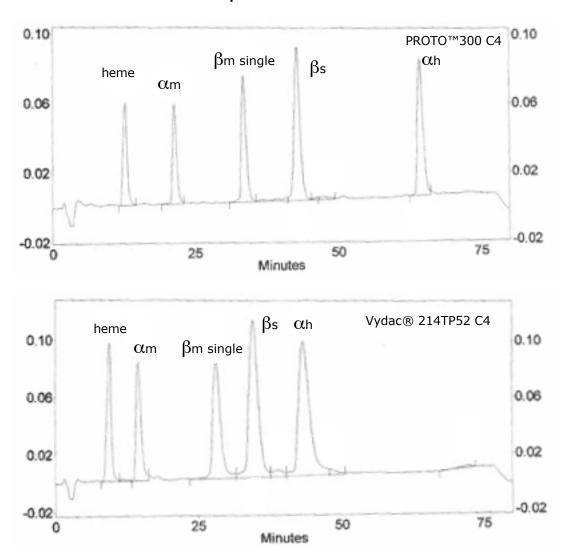


Data courtesy of Sandra M. Suzuka, Albert Einstein Hematology

Copyright © 2009 The Nest Group, Inc. Southborough MA 01772-1323

Hemoglobin Separations

Mouse Hb α , β major/minor with Human sickle Hb

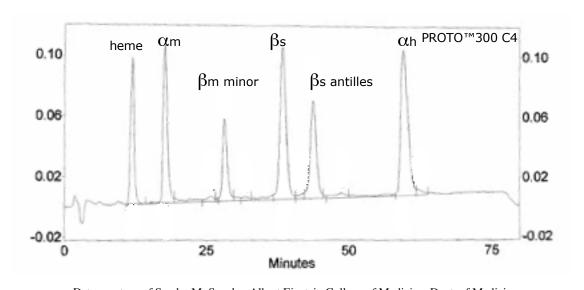

Data courtesy of Sandra M. Suzuka, Albert Einstein College of Medicine, Dept. of Medicine

Column: PROTO300 C4, H RS-2521-W045, $\ 2.1 \ x \ 250 mm, \ 5 \mu m, \ 300 \mbox{\normalfont\AA}$

Segmented Gradient: 36-56% ACN : water with 0.18% TFA

Flow Rate: 0.16ml/min

Copyright © 2009 The Nest Group, Inc. Southborough MA 01772-1323